



From Vial to Vein



# SMARTeZ™ Plus Painbloc

Elastomeric Infusion Pump

**Smart  
Pain Management**  
*Anytime • Anywhere*

**Hydrophobic Filter** in priming cap

Minimum infection risk during priming

Detachable Variable **Flow Controller**

Maximum patient safety preventing unauthorized adjustments

Ergonomic **Bolus button**

Effortless operation for both  
healthcare providers and patients of all abilities


easy to use  
no batteries  
no cords  
portable

Versatile **Applications**

- Ambulatory
- Clinical & Home settings,
- Intraoperative site
- Peripheral nerve block,
- Percutaneous administration
- Epidural route

**Epic Health EMEA**



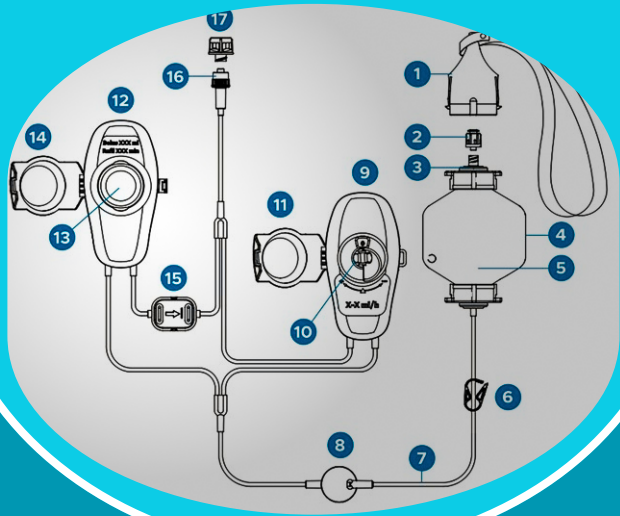

1. Hanger Clip (Only on Ref No with suffix HC)
2. Fill Port Cap
3. Top Fill Port (ISO 80369-7 Female Luer Lock)
4. Soft Cover
5. Fluid Reservoir (Elastomeric Membrane)
6. On-Off Clamp
7. Administration Tube
8. Air and Particulate Filter 1.2 micron
9. Flow Selector Module
10. Selection Knob
11. Selector Cover
12. Bolus Module
13. Bolus Button
14. Bolus Cover
15. Inline Filter 0.2 micron
16. Patient Connector / ISO 80369-7 Male Luer Lock / ISO 80369-6 NRFit (Ref No with suffix N)
17. Patient end Cap / Priming Cap for MLL / Vented Cap for NRFit

## Fixed flow

Model No (REF)	Nominal Volume	Nominal Flow Rate	Nominal Time
483020	100 ml	1 ml/h	100 hrs
483110	270 ml	5 ml/h	54 hrs
483130	400 ml	5 ml/h	80 hrs

## Fixed flow with bolus

Model No (REF)	Nominal Volume	Nominal Flow Rate	Bolus Volume	Nominal Time
485110	270 ml	2 ml/h	2 ml	15 min
485120	270 ml	2 ml/h	2 ml	30 min
485160	400 ml	5 ml/h	2 ml	15 min

## Variable flow

Model No (REF)	Nominal Volume	Variable Flow Rate
486030	270 ml	0/1/2/3/4/5/6/7 ml/h
486040	270 ml	0/2/4/6/8/10/12/14 ml/h
486050	400 ml	0/1/2/3/4/5/6/7 ml/h
486060	400 ml	0/2/4/6/8/10/12/14 ml/h

## Variable flow with bolus

Model No (REF)	Nominal Volume	Basal Flow Rate	Bolus Volume	Bolus Refill Time
487270	400 ml	0/2/4/6/8/10/12/14 ml/h	2 ml	15 min
487360	600 ml	0/2/4/6/8/10/12/14 ml/h	2 ml	15 min